
TR **TURNİKE SETİ**
Kullanma Kılavuzu

EN **TOURNIQUET SET**
Instructions For Use

REF: 330 05XX 1
330 08XX 1
330 06XX 1
330 07XX 1

CE 0044

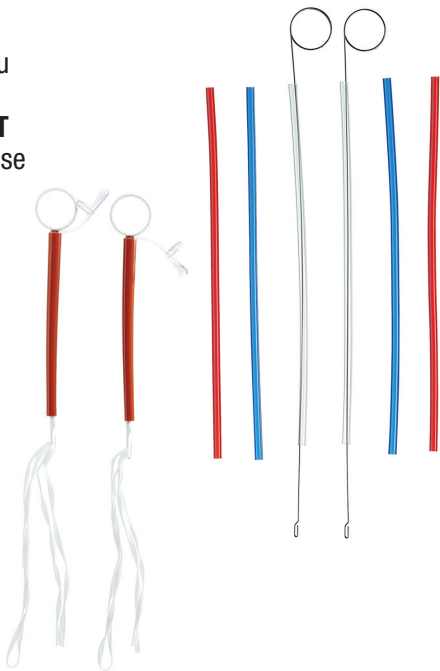SCAN FOR INFO

1. Tanım

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EN INSTRUCTIONS FOR USE

1. Description

Tourniquet Set is sterile and single use medical device that has 2 different model. CANNULATION TOURNIQUET SET consists of snare and tubing. It is used to secure a cannula. VASCULAR TOURNIQUET SET consists of snare, plug, tubing, and surgical braided tape. It is used to limit or restrict flow through a vessel.

2. Disclaimer of Liability

TOURNIQUET SET is not intended for use other than indicated in the instructions for use. The manufacturer is not responsible for problems arising from inexperience or improper use.

3. Technical Specification

The red and blue color tubes are X-ray detectable due to the radio-opaque properties of the raw material.

References	Dimension	Type of Tourniquet	Radio-Opaque Feature	Snare	Tubes (pcs x color)
330 0510 1	Tube diameter: 3.0x4.1 mm Tube length: 18 cm Snare length: 25 cm	Cannulation	Yes	2	4 x RED, 2 x BLUE, 2 x CLEAR
330 0511 1			Yes	1	1 x RED, 1 x BLUE, 3 x CLEAR
330 0512 1			Yes	2	1 x RED, 2 x BLUE, 5 x CLEAR
330 0513 1			Yes	1	1 x RED, 2 x BLUE, 3 x CLEAR
330 0514 1			Yes	1	5 x RED
330 0515 1			Yes	2	2 x RED, 2 x BLUE, 4 x CLEAR
330 0516 1			Yes	1	2 x RED, 1 x BLUE
330 0517 1			Yes	1	3 x CLEAR
330 0518 1			Yes	1	1 x RED, 1 x BLUE
330 0519 1			Yes	1	2 x RED, 2 x BLUE, 1 x CLEAR
330 0520 1			Yes	1	2 x RED, 2 x BLUE
330 0521 1			Yes	1	2 x RED, 2 x BLUE, 2 CLEAR
330 0531 1	Tube diameter: 3.0x4.1 mm Tube length: 12.5 cm Snare length: 18 cm	Cannulation	Yes	2	1 x RED, 3 x CLEAR, 2 BLUE
330 0532 1	Tube diameter: 3.0x4.1 mm and 4.0x5.0 mm Tube length: 18 cm		Yes	1	1 x RED, 2 x BLUE, 3 x CLEAR
330 0541 1	Tube diameter: 3.0x4.1 mm and 4.0x5.0 mm Tube length: 18 cm		Yes	2	2 x RED, 2 x BLUE, 4 x CLEAR
330 0542 1	Tube diameter: 4.0x5.0 mm Tube length: 18 cm Snare length: 25 cm		No	1	3 x CLEAR
330 0543 1	Tube diameter: 3.6x5.33 mm Tube length: 18 cm		No	-	4 x CLEAR
330 0544 1	Tube diameter: 3.0x4.1 mm and 4.0x5.0 mm Tube length: 18 cm		No	-	8 x CLEAR
330 0610 1	Tube diameter: 4.8x6.8 mm Tube length: 12.5 cm		No	-	5 x CLEAR
330 0710 1	Tube diameter: 3.0x4.1 mm Tube length: 18 cm Snare length: 25 cm		Yes	2	2 x RED, 2 x BLUE, 4 x CLEAR
330 0590 1	Tube length: 13 cm	Vascular	No	-	2 TUBES
330 0591 1			No	1	1 TUBE
330 0592 1			No	1	5 TUBES
330 0593 1			No	-	1 TUBE
330 0594 1			No	-	5 TUBES
330 0841 1	Tube diameter: 4.8x6.8 mm Tube length: 15 cm	Cannulation	No	-	5 x CLEAR

4. Mode of Action

TOURNIQUET SET controls blood flow by applying circumferential pressure, compressing the underlying tissues and blood vessels. This pressure occludes arteries and veins, effectively stopping blood flow to the area distal to the tourniquet.

5. Intended Use

VASCULAR TOURNIQUET SET is used in a surgical procedure to secure a cannula or to limit or restrict flow through a vessel. CANNULATION TOURNIQUET SET is used to secure a cannula within such vessel in cardiopulmonary bypass surgeries.

6. Intended User

The product is designed to be used by specialist doctors. Proper surgical procedure and technique is the responsibility of the attending medical professional (physician). This guideline is furnished for informational purposes only. Each doctor must evaluate the appropriate use of this device, case by case, based upon medical training, experience and the type of surgical procedure employed.

7. Intended Environment

It is used in direct contact with tissue and blood.

8. Intended Patient Population

All ages.

9. Information on Accessories and Other Devices

Tourniquet Sets can be used with surgical cannulas.

The Cannulation Tourniquet Set is used with surgical sutures.

10. Material

CANNULATION TOURNIQUET SET: PVC, Stainless Steel.

VASCULAR TOURNIQUET SET: PVC, Cotton, PP.

11. Indications

TOURNIQUET SET is indicated in the surgical procedures where it is desired to limit or restrict flow through a vessel.

12. Contraindications

TOURNIQUET SET should not be used in patients with known hypersensitivity to PVC, PP, Cotton, stainless steel components.

13. Residual Risk

According to the risk analysis carried out by the Bıçakçılar Risk Management Team, the residual risk is determined to be as far as reasonably possible.

14. Recommendations For Use

- Remove the Tourniquet Set from the packaging with aseptic technique.
- To shorten the tubing, remove the snare and aseptically crop the tourniquet tubing to desired length. Reinsert snare into the tubing.

CANNULATION TOURNIQUET:

- Place the surgical sutures around the cannulation area.
- and hook around the wire snare.
- Pull snare through the tube.
- Check if the required tension is obtained and clamp the suture end with a hemostat.
- To remove the cannula applied on a vessel, cut the sutures around the vessel, open the clamp and remove the sutures while pulling the tube.
- Discard the tourniquet and the tube properly.

VASCULAR TOURNIQUET:

- Position tourniquet around cannulated vessel and hook tourniquet with snare.
- Pull snare and tourniquet through tube and unhook tourniquet from snare.
- While applying tension to tourniquet, with one hand, slide tubing down the tourniquet. Once desired tension has been applied, push retaining plug into tube.
- Make sure that retaining plug is well seated and tourniquet is applying desired tension.
- Wrap excess tourniquet material around end of tube and cinch through slot in retaining plug to hold material securely.
- To remove tourniquet, release and unwrap excess material, unplug tube and slide off tourniquet.
- Discard the tourniquet and the tube properly.

DURING THE APPLICATION

- Apply tension gradually to the loop to achieve controlled occlusion of the vessel.
- Use the minimum tension necessary to achieve effective hemostasis.
- Secure the TOURNIQUET in position using the fixation mechanism to maintain stable occlusion.
- Continuously verify the effectiveness of the occlusion by assessing bleeding control.
- Monitor the condition of the vessel and surrounding tissues throughout the occlusion period."
- Assess distal perfusion where clinically relevant."
- Regularly evaluate the duration of occlusion to avoid prolonged ischemia.

RELEASE OF OCCLUSION

- Release the occlusion in a controlled and gradual manner, unless rapid release is clinically required.
- Observe the vessel and surgical field during reperfusion for any abnormal bleeding or reactions.

15. Warnings

- Inspect the product for damage prior to use. Damage to the product (i.e., kink crimps, distortions or cuts) may result in obstruction or disruption of flow during use and product performance may be affected. Do not use damaged product.
- If the packaging is damaged, do not use as sterility of the device may have been compromised.
- Use aseptic techniques during the procedure.
- A spare product must always be available during operation.
- Make sure that surgical tape is at appropriate stress. Excess force may tear the surgical tape apart.
- Excessive force during the application may lead to deformation on the wire or separation of the stylet tip. Check the product integrity during the procedure.
- The product is produced for single use and for single patient use only: the product comes in contact with human blood, body fluids, liquids or gases.
- Re-sterilization of the product may negatively affect the product performance. Therefore, reuse on other patients may cause cross-contamination, infection and sepsis. Do not resterilize.
- Do not apply excessive traction, as this may result in vascular injury or damage to surrounding tissues.
- Do not exceed the recommended occlusion duration, as prolonged occlusion may lead to ischemia and tissue damage.
- Apply occlusion pressure gradually and only to the extent necessary to achieve effective hemostasis.
- Use caution in calcified, fragile, or diseased vessels, as effective occlusion may not be achieved.
- Rapid release of occlusion may cause sudden reperfusion and associated physiological effects.
- Monitor the vessel and surrounding tissues continuously during occlusion.
- Ensure correct positioning of the tourniquet around the vessel before applying tension.
- Do not use the device on vessels with severely compromised integrity, including fragile, dissected, or extensively diseased vessels, unless clinically justified.

16. Adverse Events and Complications

Hemostasis, Vascular wall damage, and hemodynamic changes.

17. Undesirable Side Effects

Local vessel injury.

18. Caution

Use of tourniquets in neuropathies, fragile vascular structures, systemic lupus erythematosus, underlying coagulation disorders, and cachectic patients must be done with caution.

19. Packaging and Transportation

The primary package consists of snare, plug, tubing, and surgical braided tape. It comprises a ready-to-use sterile packaged Tourniquet Set.

20. Expected Clinical Benefits

The use of the TOURNIQUET SET helps to secure the cannulas, to sign the veins and to create a bloodless field in surgery.

21. Performance Requirements and Limitations

VASCULAR TOURNIQUET SET must be monitored to ensure it is used for the shortest time necessary. It is not recommended to use a tourniquet for more than 2 hours since after this time of ischemia, progressive venous acidosis will occur in the region the tourniquet effects.

CANNULATION TOURNIQUET SET can be used for 6 hours during cardiopulmonary bypass surgery.

22. Disposal of The Device

The product should be discarded as contaminated medical waste.

23. Notice Regarding Serious Incident

Any serious incident that will occur in relation to the device should be reported to Biçakçılar and the correspondent competent authority of the Member State in which the user and/or patient is established.

24. Basic UDI-DI

8695476DP330050RY

25. Shelf Life

The shelf life of the device is 5 years.

26. Storage

The device is recommended to be kept in its original packaging and stored at a temperature between 5°C - 35°C until use. Protect from direct exposure to sunlight. Under proper storage conditions, the device may be used up to the expiry date indicated on the label.

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